


## TECHNICAL INFORMATION

Catalog Number: 150109, 191406, 194623, 194806

### Ammonium Chloride

**Structure:**

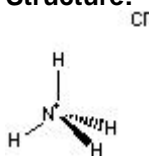

**Molecular Formula:** NH<sub>4</sub>Cl

**Molecular Weight:** 53.5

**CAS #** 12125-02-9

**Synonyms:** Ammonium chloride; Amchlor; Ammoneric; Darammon; Sal ammoniac; Ammonium muriate; Salmiac

**Physical Appearance:** White crystalline powder

**Solubility:** Soluble in water (100 mg/ml - clear, colorless solution). The presence of hydrochloric acid or sodium chloride decrease water solubility.<sup>1</sup> Also soluble in methanol and ethanol.<sup>1</sup>

**pH:**<sup>1</sup>

1% aqueous solution: 5.5

3% aqueous solution: 5.1

10% aqueous solution: 5.0

**Description:** Ammonium chloride is used in both industrial and research applications. Industrial applications include as a flux for coating sheet iron with zinc; tinning; in dry and Leclanch batteries; dyeing, freezing mixtures, electroplating, to clean soldering irons, safety explosives, lustering cotton, tanning; in washing powders; manufacturing of dyes; in cement for iron pipes; for snow treatment (slows melting on ski slopes).<sup>1</sup>

In research, ammonium chloride has been used in:

- To lyse human red blood cells<sup>2,5,8</sup>
- Used in the study of basic calcium phosphate crystals in fibroblasts<sup>4</sup>
- Used in the isolation of proteins from 50S ribosomal subunits of *Bacillus stearothermophilus*<sup>3</sup>
- Used in biological buffers

**Availability:**

Catalog Number	Description	Size
150109	Ammonium Chloride	100 g 500 g 1 kg
194623	Ammonium Chloride, cell culture reagent	100 g 500 g 1 kg
194806	Ammonium Chloride, molecular biology reagent	100 g 500 g 1 kg
191406	Ammonium Chloride ACS Reagent Grade	500 g 1 kg 5 kg

**Also Available:**

Catalog Number	Description	Size
----------------	-------------	------

540113	<a href="#">Ammonium Chloride-<sup>15</sup>N, purity 99% <sup>15</sup>N atom</a>	1 g
--------	--	-----

#### References:

- *Merck Index*, **12th Ed.**, No. 537.
- Baerlocher, G.M., et al., "Telomere length measurement by fluorescence in situ hybridization and flow cytometry: tips and pitfalls." *Cytometry*, **v. 47(2)**, 89-99 (2002).
- Gewitz, H.S., et al., "Reconstitution and crystallization experiments with isolated split proteins from *Bacillus stearothermophilus* ribosomes." *Biochem. Int.*, **v. 15(5)**, 887-895 (1987).
- Halverson, P.B., et al., "Intracellular calcium responses to basic calcium phosphate crystals in fibroblasts." *Osteoarthritis Cartilage*, **v. 6(5)**, 324-329 (1998).
- Kang, E.M., et al., "Mobilization, collection, and processing of peripheral blood stem cells in individuals with sickle cell trait." *Blood*, **v. 99(3)**, 850-855 (2002).
- Korber, P., et al., "HSP15: a ribosome-associated heat shock protein." *EMBO J.*, **v. 19(4)**, 741-748 (2000).
- Locatelli, C., "Measurement of voltammetric peak area and resolution of overlapping peaks in the simultaneous determination of copper, lead, cadmium, and nickel in environmental matrixes." *J. AOAC Int.*, **v. 83(6)**, 1321-1326 (2000).
- Ridings, J., et al., "Purification of cord blood lymphocytes." *J. Immunol. Methods*, **v. 195(1-2)**, 43-48 (1996).